

m y s é

English User's Manual come with this product.

Please download the other language's User's Manual version from this link (<https://www.ya-man.co.jp/en/support/documents/>) or scan with QR code with your camera phone.

Please carefully read User's Manual to ensure safe and proper use of the device.



TO ACCESS SCAN WITH
YOUR CAMERA PHONE

User's Manual

Home-care Beauty Device

m ys é Deep Skin Clear

MS-43

Introduction

Thank you for purchasing this YA-MAN product.

- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Be sure to include this User's Manual if you are transferring ownership of the device to another person.

How to Use

1

Turn on the device

Hold down the POWER/MODE button
(for approx. 3 seconds).

PEELING MOIST BODY
C C C



2

Select the mode

Press the POWER/MODE
button.

PEELING MOIST BODY
C C C



3

Select the level

Press the LEVEL button.

PEELING MOIST BODY
C C C



4

Place the device on your skin

Use the device on the wet skin.

Table of Contents

Safety Precautions	p3
Names of Components	p10
Before Use	p12
How to Use	p18
After Treatment	p24

Additional Information	p27
After-sales service	p31
● Product Specifications	p32
● Warranty	p34

■ Safety Precautions

- Read these Safety Precautions carefully before using the device.
- Use the device sensibly in accordance with the prescribed steps.
- Immediately discontinue use and arrange to have the device repaired if any abnormality occurs.
- Do not use the device in a defective state, as doing so could cause a physical harm to the user.
- The precautions listed in the User's Manual are very important. For your safety, please follow them carefully.
- Precautions are classified as DANGER, WARNING, or CAUTION in order to clearly specify the level of danger and possible injury or damage that may occur if the device is misused.

Classifications of Precautions



Improper use of this device poses a severe and likely risk of death or serious injury.



Failure to operate the device according to this type of precaution could result in death or serious injury.



Improper use of this device may result in an injury or a property damage*.

Symbols



This symbol shows you what you must not do.



This symbol shows you what you must do.



This symbol shows you what you must be careful with.

*'Property damage' refers to damage affecting dwellings, household belongings, or domestic animals and pets.

DANGER (Use and storage)

- The following individuals should not use the device:

Failure to comply may cause electronic medical devices to malfunction.

- ◆ Do not use the device in conjunction with any medical electronic equipment, including:
 - pacemakers or other implanted electronic medical devices;
 - artificial heart-lung machines or other electronic life-support medical devices; and
 - electrocardiographs or other worn electronic medical devices.



Prohibited

WARNING (Use and storage)

- Do not use under the following conditions:

Failure to comply may cause accidents or other problems.

- If you are under a doctor's care · If you are in a physical growth stage · If you have an acute disease · If you have tuberculosis · If you have a fever
- If you are suffering from hypertension · If you have an infectious disease · If you have a malignant tumor · If you have a blood disease
- If you have a visceral disease · If you have a skin disease · If you have an allergic constitution · If you have atopic dermatitis · If you have skin problems
- If you are severely fatigued · If you have sensitive skin · If you are physically unwell · If you are under the influence of alcohol
- If you are pregnant, menstruating, or lactating · If you take medication that induces drowsiness
- If you have unhealthy skin blemishes caused by physical irritation · If you are experiencing skin irritation or redness due to sunburn
- If you have capillary dilatation due to long-term use of steroid hormonal drugs or a liver disorder · If you have excessive myopia (-6D or higher)
- If you had an eye or face operation; · If you have diseases of the eyes or face or you are under treatment



Prohibited

- Do not use the device on the areas indicated below or outside of specified areas.

Failure to comply may cause an accident or other problems.

- Cuts · Areas where surgery (including plastic surgery) was performed · Areas where cosmetic treatment was performed
 - Anywhere inside the mouth, around the pubic region, or any mucous membrane · Acne · Areas where skin suffers from purulent inflammation
 - Areas with itchiness or hot sensations · Bulging moles · Viral warts · The Adam's apple · Eyeballs · The eyelids or around the eyes
- If you are uncertain as to whether it is safe to use the device, consult a physician prior to use.

WARNING (Power supply)

- Do not wash or submerge the charging base cord in water.

Failure to comply may cause an electric shock, short circuit or fire.



Exposure
to water
is prohibited

- When the device is not in use, make sure to turn off the power, hold the head of cord and disconnect the AC adapter from the power outlet. Do not insert or disconnect the plug with wet hands.

Failure to comply may cause an electric shock, short circuit or fire.

- In the unlikely event that an abnormality arises, turn off the power.

Do not use the device if smoke, abnormal noises, or an odor is detected, or if it is otherwise operating in an abnormal state.



Mandatory
compliance

- Make sure that the charging base cord is fully inserted into the USB port.

Failure to comply may cause a fire, electric shock, or failure of the device.



Mandatory
compliance

- Do not use the device if the charging base cord is damaged or if the USB port is loose.

Failure to comply may cause an electric shock, short circuit or fire. Contact your local distributor for advice on repairing the device.



Use of a
damaged cable
or plug is
prohibited

- Be sure to use the dedicated charging base for this device. Plug the AC adapter (commercially available) directly into a wall outlet and do not overload the electric circuit.

Failure to comply may cause an electric shock, short circuit or fire.



Overloading
of an electric
circuit is
prohibited

- If any problem occurs during use, including the presence of an abnormal odor or excess heat, stop using the device.

Failure to comply may cause the device to become damaged or catch a fire.



Do not use
when the
device is hot

- Do not place the device near fire or flames, and do not dispose of the device in a fire. Do not use, charge, or store in a high-temperature area or if the device is hot. The device features a built-in rechargeable battery. Failure to comply may lead to the generation of heat, a fire or explosion.



Prohibited

WARNING (Use and storage)

- Do not use the device for any purpose other than those stated in the User's Manual.

Failure to comply may result in an injury or malfunction of the device.



Non-intended uses are prohibited

- Do not damage, break, modify, bend by force, pull, twist, or bundle the charging base cord.

Do not place heavy objects on the charging base cord or allow it to get pinched in between other objects because this may damage the cable and cause a fire or electrical shock hazard.



Careless or reckless use is prohibited

- Do not place the device on an unstable location. Do not drop or strike the device.

Doing so may cause a failure of the device. An injury may also be caused if the device is dropped on a foot or other exposed body parts.



Dropping of the device is prohibited

- Do not use or leave the device near fire, in direct sunlight, or in a dusty location. Do not bring the device close to a heater.

Failure to comply may lead to the generation of heat, a fire or explosion.



Fire hazard

- Do not disassemble, attempt to repair or modify this device in any way.

Failure to comply may cause a fire, electric shock, or an injury. Contact your local distributor for advice on repairing the device.



Disassembly is prohibited

- Do not drop the device into water.

Do not immerse the device in water. Do not use or store the charging base in the bathroom or other places where there is high humidity.

Failure to comply may cause an electric shock, short circuit or fire.



Prohibited

- Do not allow persons who are unable to express their intentions on their own to use the device. Use by unattended people with physical disabilities should be avoided.

- Keep away from children and do not allow them to use the device.

Failure to comply may result in an accident or injury.



User restrictions

CAUTION (Use and storage)

- If you are currently undergoing dental treatment or have teeth fillings, choose a lower intensity level and apply the device on a different area.

You may experience unpleasant sensations in your teeth due to the ion/ultrasonic stimulation. Choose a lower ion/ultrasonic level and apply the device on a different area.



Mandatory compliance

- Discontinue use if you start to feel fatigued or you see any skin reaction during use.

If symptoms persist the following day and beyond, consult a physician.



Mandatory compliance

- Be sure to use the device only for the length of time prescribed.

Excessive use may cause skin damage.



Mandatory compliance

- If you feel an unpleasant tingling on your skin, move to another area. You may also reapply water or lotion.

Failure to comply may cause skin damage.



Mandatory compliance

- When using the device in treatments, move it continuously over the skin and do not stop over any one area.

Stopping the device over one area during treatment may cause skin damage.



Mandatory compliance

- Do not press too hard against your skin.

Failure to comply may cause skin damage.



Mandatory compliance

- Reapply water or lotion during treatment to ensure that there is always some water or lotion between the treatment head and the skin.

Failure to comply may cause skin damage.



Mandatory compliance

- After each use, refer to 'After treatment' (page 24) and wipe off any dirt or moisture that have adhered to the device.



Mandatory compliance

- Remove any accessories (necklace, belt, rings, earrings, etc.) before using the device.

Failure to comply may result in an accident or injury.



Prohibited

CAUTION (Use and storage)

- Do not use the device excessively on the same skin area. Excessive peeling may cause skin damage.



Prohibited

- Be sure to perform a pretest before using the device for the first time. Use only level LOW upon first use until you are accustomed to use.



Prohibited

Failure to comply may result in an accident or injury. Be sure to adjust the level only after your body has become accustomed to the device.

- Do not allow the device to come into contact with another beauty device or electronic appliance while the power is on.



Prohibited

Failure to comply may result in an injury or malfunction of the device.

- Do not use the device on pets. Do not use the device when pets are around or close to you.



Prohibited

Pets may react to the operating sound of the ultrasonic waves and may panic.

- Before using the device, check that the head has no scratches, dirt, or other abnormalities. Do not use the device if the head has any scratches or damage.



Prohibited

Failure to comply may cause skin problems or a malfunctioning of the device.

- Do not use with any substance other than cosmetic products (except those not appropriate for the use) and water.

Failure to comply may cause skin problems or a malfunctioning of the device.

<Cosmetics not appropriate for use>

- Oil
- Scrub
- Stimulating toner
- Facial cleanser



Prohibited

- Do not use the device with water or lotion on your hands. Your hands may slip, causing unforeseen accidents, or the gel/water may enter the spaces between the main unit and buttons, causing the device to malfunction.



Prohibited

- Push the head firmly all the way into the back of the main unit to fix it in the place.

If the head is not fixed firmly in place, the ultrasonic wave output may weaken or the device may malfunction.



Prohibited

- Do not wipe the device with alcohol, thinner, benzene, nail-polish remover, or any other solvents.

Failure to comply may cause damage to the surface of the device.



Prohibited

CAUTION

- Sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children shall not play with the device.
- Cleaning and user maintenance shall not be made by children without supervision.
- Cleaning and user maintenance shall not be made by children without mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.



Mandatory compliance

WARNING (when using outside of Japan)

- Use only with a 5 V DC USB power supply (warranty only available in Japan).
Failure to comply may cause a fire or electric shock.



Mandatory compliance

- Be sure to check the plug configuration and the voltage of the power outlet before use.
Failure to comply may cause the generation of heat, a fire, or a malfunction of the device.




Fire hazard

- The plug configuration and voltage may vary depending on the region and the type of accommodation.
Always be sure to check for compatibility before using the device.
- For information about plug adapters, consult a consumer electronics retailer.


Names of Components

Main Body

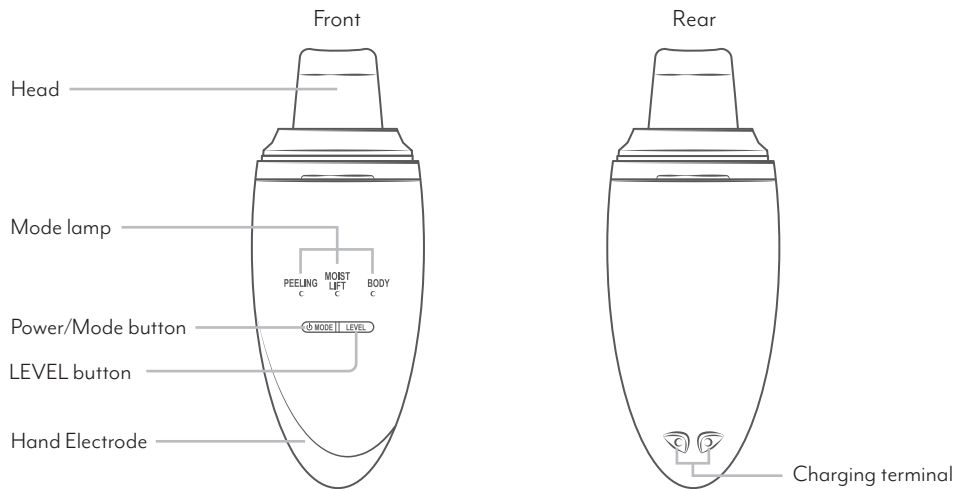

* Illustration is when the wide head is mounted

Accessories

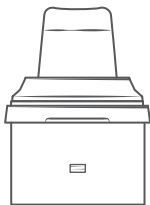


Wide head*

* Attached to the main unit.




Spoon head

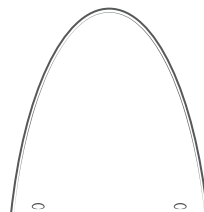

Cover common for Wide head
and Spoon head: One

User's Manual

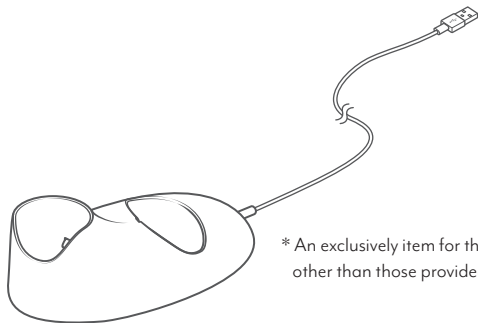

Charging base (with USB cable)


* An exclusively item for this product. Do not use attachments
other than those provided with the product.

* The design and specifications of the device and accessories are subject to change for improvements without prior notice.


■ Before Use

Charge the device

Using a household power outlet

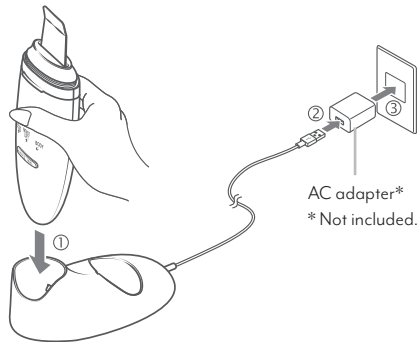

Using a computer

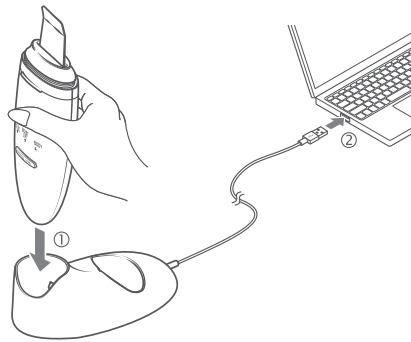

- Place the device so that it touches the charging terminal of the charging base.
- Fully charge the battery before using the device for the first time after purchasing. When using the device for the first time or when the device has not been used for an extended period of time, it is possible that you will be unable to use the device for very long even with a full charge displayed. In such a case, recharge the device before using it.

Charging	The Mode lamps light up in sequence.
When charging is complete	All of the Mode lamps light up.
Insufficient charging	All the lamps flash. Charge the battery.
Charging time	Approx. 1.5 to 2.5 hours (at room temperature 28°C). The actual charging time will differ depending on the remaining battery level, the charging environment, and the AC adapter (commercially available) being used.
Operating time	Approx. 30 minutes. The actual operating time will differ depending on the remaining power level and the charging environment.

CAUTION

- The charger base is not waterproof. If there are any water droplets on the charging terminal, make sure to wipe dry thoroughly before use. Charging while water droplets are remaining on the device could cause a failure.

Functions

Auto Power-off

The device will switch off automatically around 5 minutes after a mode starts to operate (switching the mode will reset the Power-off timer, but changing the level in the same mode will not reset it).

Water resistance class

This product is IPX5 level waterproof (JIS waterproof protection class).

* The device is not completely waterproof. Do not submerge the device in water or leave it continuously under running water for 3 minutes or longer.

Areas of use

This product comes with two types of heads. Use different heads for different areas of your body.

Wide head: Face, elbows, knees, and heels.

Spoon head: Face (recommended for nose, between the eyebrows and around the chin.)

Switching modes

The device has three different treatment modes. It is possible to switch modes as follows by pressing the POWER/MODE button with the power on status. If you want to turn the power OFF during use, hold down the POWER/MODE button (approx. 3 seconds).

- Double Peeling Mode (Ultrasonic + Ion Cleansing)
Both Wide head and Spoon head can be used.
- Moist Lift Mode (Ultrasonic + Ion Moisture)
Use the wide head
- Body Mode (Ultrasonic)
Use the wide head.

* If the head is different, it will not be possible to switch modes.

Adjusting the level


The device has two operating levels, LOW and HIGH.

Each time the Level button is pressed, the level is changed.

Be sure to start using the device with level LOW.

* The ultrasonic waves become stronger when the level is increased.

* The ion level for each mode cannot be changed.

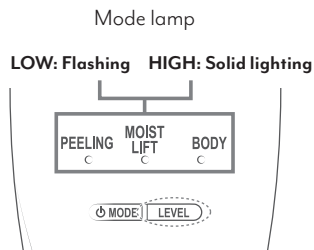

Pretesting on the skin

When using the device or a skincare product for the first time, always test for compatibility with your skin to prevent skin problems.

1 Clean your skin.

Wash the inner side of your arm with a soap and water, then wipe dry with a towel.

2 Apply water or lotion to the inside of your arm and use the device while moving it for about 1 minute in Double Peeling Mode.

- Use level LOW.
- Do not rub the pretested area.
- Use the wide head.

3 Continue to use the product only if you experience no reddening or burning sensation in your skin within 24 hours.


CAUTION

- If a rash, reddening, itchiness, stinging sensation or other abnormalities occur at the pre-testing area, immediately stop the pretest and wash the area.
- Check carefully to make sure that there are no skin abnormalities before using the device. Stop using the device if the treatment does not agree with your skin.

How to Detach/Attach the Head

Attach either the Wide head or the Spoon head before use depending on the area of use.

Please see pages 18 to 23 for skin areas that can be treated using the device.

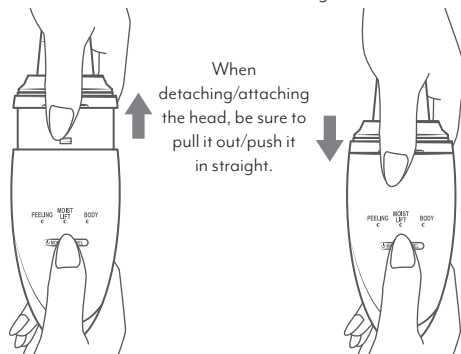

- Thoroughly wipe off and remove any moisture on the main unit and head.
- Align the orientation of the main unit and the head and push the head **STRAIGHT** into the main unit as far as it will go. Both fittings have to face each other properly. Otherwise, the attachment cannot be installed. Check the direction in which the head should be inserted before attaching it.
- When detaching the head, grab the head and pull out straight. If the head is loose, the power supply may cut off partway through use. If you pull out in the oblique direction, the fitting parts may be deformed or damaged.
- Check whether there are scratches or dirt on the head or connecting part. Failure to comply may cause skin problems or a malfunctioning of the device.

CAUTION

■ How to Use

Recommended use

Use the device about one or two times per week.

As a general rule, use the device for 5 minutes or less each time for one area in one mode (10 minutes for Moist Lift Mode).

* This is only a guideline. Please observe your skin condition and use the device accordingly. Excessive use may cause skin damage.

* The device will switch off automatically approximately 5 minutes after a mode starts to operate (Auto Power-off function). If you want to continue using the device, turn it on again.

Preparation

1 Clean your skin.

- Remove oil and dirt from your skin so that it is in a clean condition prior to the treatment.
- Make sure to remove makeup before using the device on your face.

3 Check that your skin has no abnormalities, such as pimples or breakouts.

2 Select the head and attach it to the main unit.

- Do not use the device on skin areas where there are abnormalities.
- Do not use the device on the areas around the eyes or eyelids.
- If you perceive any skin abnormalities during use, immediately stop using the device and consult a dermatologist.

How To Use The Device In Each Mode

● Common to all modes

1 Thoroughly apply water or lotion to your skin with an atomizer.

* Apply lotion when using Moist Lift Mode.

2 Hold down the **POWER/MODE** button (approx. 3 seconds) to turn on the device.

Press the button again briefly for Moist Mode or press it briefly twice for Body Mode.

3 Select your preferred level.

4 Treatment.

-
- Reapply water or lotion during treatment to ensure that there is always some moisture between the treatment head and the skin.

Using the device on dry skin may cause skin damage.

- Move the head continuously over the skin and do not stop over any one area.

Stopping the device over one area during treatment may cause skin damage.

- Do not use with any substance other than lotion (except those not appropriate for the use) and water.

Failure to comply may cause skin problems or a malfunction of the device.

<Cosmetics not appropriate for use> · oil · scrub · stimulus toner · facial cleanser

CAUTION

Points to note when using the device and how to place the device on the skin

- After thoroughly applying water or lotion to the face using an atomizer, place the device gently on the skin and slowly slide it over the skin to achieve effective results. If there is insufficient water on the skin, reapply the water to the skin.
- Use the device after warming the skin with a hot steamed towel or after bathing when the skin is still warm to achieve effective results.
 - * Do not use the Body Mode on the face.
 - * When the Spoon head is attached, Moist Lift Mode/Body Mode will not start.

CAUTION

- Do not drop the device into water.
- Do not use the device with water or lotion on your hands.

Failure to comply may result in a malfunction of the device, an accident or injury.

● Double Peeling Mode


Place the device gently on the skin then move it slowly in the direction of the arrows.

● After using the device, apply lotion or serum to the skin.

* Move the device continuously over the skin and do not stop over any one area.


How to use the Spoon Head

· Moving in circles

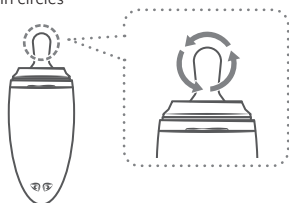

Apply the device so that the whole area shown by the dotted line contacts your skin and move the device in circles.


· Moving backward and forward

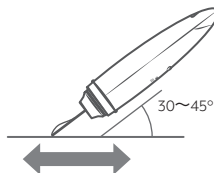

Recommended areas of use:
nose, between the eyebrows,
around the chin




How to use the Wide Head

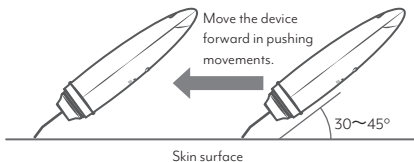


Areas you should NOT treat

Direction in which to place the device

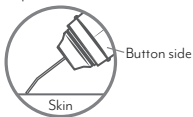

Move the device towards the center of the face, as going against the pores.

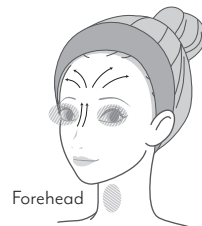


Move the device in an upwards direction, covering the entire face.

● Moist Lift Mode (with Wide Head Only)

As shown in the diagram, place the device gently on the skin at an angle of approx. 20° to 30° , then move it slowly in the direction of the arrows.


- * Check the front and back of the head before placing it on your skin.
- * Move the device continuously over the skin and do not stop at any one area.
- * When Spoon head is attached, the Moist Lift Mode will not start.

How to use

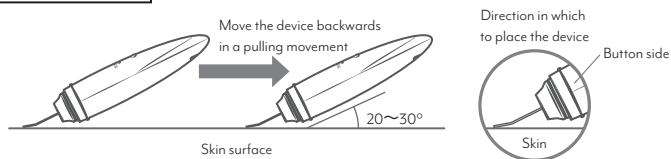

Move the device upwards along the face contour

 Areas you should NOT treat

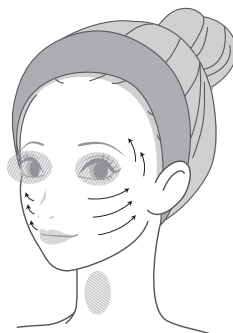


● Body Mode (Wide Head Only)


As shown in the diagram, place the device gently on the skin at an angle of approx. 30 to 45°, then move it slowly in the direction of the arrows.

- After using the device, apply lotion or serum to the skin.
- * Check the front and back of the head before placing it on your skin.
- * Move the device continuously over the skin and do not stop at any one area.
- * When Spoon head is attached, the Body Mode will not start.

How to use

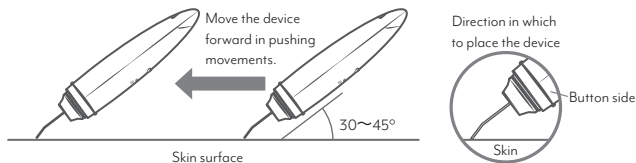


● Elbows

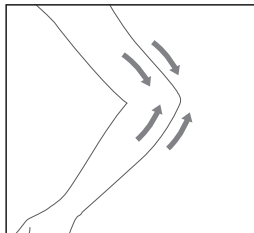


Bend the elbow slightly and move the device from the outside to inwards, covering the entire elbow.

● Knees

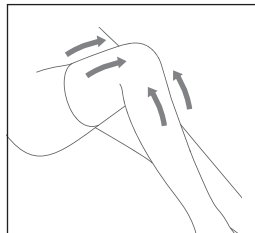

Bend the knee slightly and move the device from the outside to inwards, covering the entire knee.

● Heels

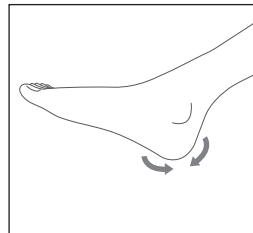

Move the device from the outside to inwards, covering the entire heel.

■ After Treatment

How To Clean



Clean the device after each use for hygienic reasons.

1 Make sure the device is turned off.

2 Detach the head from the device.

- For how to detach the head, refer to 'How to Detach/Attach the Head' (page 17).


3 Wash off any grime on the head under running water

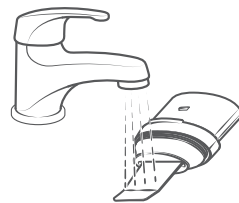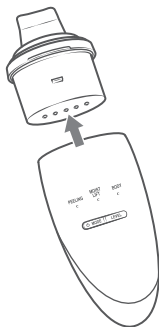

CAUTION

- When cleaning the device, be sure to turn the power OFF.
Failure to comply may cause an electric shock, short circuit or fire.

- 4** Wipe the surface of device and connecting terminals of the treatment head and main unit with a dry cloth, such as a towel or cotton swab to remove all grime and moisture.

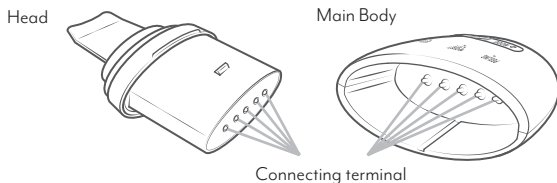

- 5** Dry the device in a location with good ventilation.

CAUTION

- Do not wipe the device with alcohol, thinner, benzene, nail-polish remover or any other solvents.
- NEVER leave the main unit or the head in water. Doing so may cause a failure of the device.
- Be sure to wipe off moisture after use. Also, do not leave the device in hot or humid places, such as in a bathroom.

Leaving the main unit wet may cause the charging terminal to corrode.



Storing

- When not using the device for an extended period of time, be sure to carry out maintenance on the main unit and ensure that the cover is on the head before storing the device.
- Be sure to thoroughly remove any moisture on the main unit or head before storing the device.
- When not in use, turn off the power, unplug the AC adapter (commercially available) from the power outlet, remove any moisture from the device, and dry well before storing.
- Do not leave the device near fire or flame, under direct sunlight, or in humid places such as the bathroom. Store only in a clean, well-ventilated places.
- Keep the device out of the reach of children.

Maintenance and Inspections

- If you will not use the device for an extended period, be sure to inspect it carefully and check that the device works properly and safely before restarting use.
- If the device is malfunctioning, go to 'Checklist for Possible Problems' (page 27 to 28) and attempt to identify the problem before arranging to have the device repaired. If the problem persists, contact your local distributor, making sure that you have your purchase details ready when you call.
- If the device malfunctions, appears damaged or needs service, please contact your local distributor.

■ Additional Information

Checklist for Possible Problems

If you perceive any abnormality with the device during use, check to see whether any of the following symptoms are present.

Symptom	Possible Causes and Solutions
The device does not turn on	● Hold down the POWER/MODE button (approx. 3 seconds) to turn on the power.
	● There is a build-up of grime or water on the connecting parts of the treatment head and/or main unit. Remove grime and water using a towel, cloth or cotton bud (page 24).
	● Charge the battery (page 12).
The device only operates for a few minutes even after charging	● Charge the battery (page 12).
	● There is a build-up of grime or water on the connecting parts of the treatment head and/or main unit. Remove grime and water using a towel, cloth or cotton bud (page 24).
	● The battery has expired. Inquiry the Global After Service (https://www.ya-man.co.jp/en/support/repair/) or your local distributor.
The device does not start charging	● Insert the AC adapter (commercially available) or charging base cord properly (page 12).
	● The main unit is not set firmly on the charging base. Insert the main unit into the charging base correctly (page 12).

Symptom**Possible Causes and Solutions****Ultrasonic waves are painful**

- The head is touching the skin in a specific area. Please apply the head correctly (pages 22 to 23).
- Lower the level using the Level button.
- Device is rubbed too strongly against the skin. Slide the device lightly over the skin.

Ultrasonic waves are weak

- The head is not touching the skin correctly. Please apply the head correctly (pages 22 to 23).
- The head is not attached correctly. Be sure that the head is firmly inserted all the way into the main unit (page 17).

The device shuts off

- The Auto Power-off function (approx. 5 minutes) worked. This is not a malfunction. If using the device continuously, turn on the power again (page 14).
- The head is not attached properly. Be sure that the head is firmly inserted all the way into the main unit (page 17).
- There is a build-up of grime on the connecting parts of the treatment head and/or main unit. Remove grime and water with a towel, cloth or cotton bud (page 24).

All of the lights are flashing

- Weak battery. Charge the battery (page 12).

Frequently Asked Questions

Q

Why do I need to moisturize the skin area being peeled?

A

The skin area being peeled needs to be moisturized so that the ultrasonic waves are carried evenly to the skin. **NEVER** use the device on dry skin.

Q

What is the recommended duration and frequency of use?

A

- Use the device about 1 to 2 times a week.
As a general rule, use the device for 5 minutes or less each time for one area in one mode (10 minutes for Moist Lift Mode).
* This is only a guideline. Use the device according to your skin condition.
* The device will switch off automatically approximately 5 minutes after a mode starts to operate (Auto Power-off function). If you want to continue using the device, turn it on again.
- Check carefully to make sure that there are no skin abnormalities before using the device. Stop using the device if the treatment does not agree with your skin.

Q

My skin becomes red

A

Rubbing the head strongly against the skin may cause skin damage. Be sure to slide the head lightly over the skin. Also, reapply water or lotion during treatment to ensure there is always gel or water between the treatment head and the skin.

Q

The main unit is difficult to move

A

If the skin surface is dry, the head may not be able to move over the skin smoothly. For effective treatment, apply water or lotion to your skin before using the device.

Q Can I use the device in the bathroom?

A The device may be used in the bathroom, but please note that it cannot be held under running water for long periods of time or submerged in water.

Q Will I feel a tingling sensation on the skin?

A Because of differences between individuals, some people may feel a stinging sensation from the weak electric current. If you feel any skin discomfort, stop using the device and consult a dermatologist. Note that if the treatment head is not touching the skin properly, the ultrasonic waves may cause a tingling sensation.

Q Can I use the device over makeup?

A The device is not a makeup remover, so be sure to remove all makeup before using the device.

Q When using the device, can I only apply water to the skin?

A Lotion may also be used with the device, but do not use scrubs or other products that irritate the skin.

However, do not use any of the following as they may cause skin problems or a malfunction of the device:

<Cosmetics not appropriate for use>

- Oil
- Scrub
- facial cleanser
- Stimulating skin lotions and other cosmetics

Q There is no atomization or spray

A Check whether the head is firmly inserted all the way into the main unit. Detach the head and then reinsert it.

Depending on the water volume, angle, and type of lotion, the water may not spray efficiently. If the plates are grimy (especially with oil), the water may not spray efficiently. Thoroughly remove any grime on the plates.

Depending on the water volume, angle, and type of lotion, the water may not spray efficiently.

■ After-sales service

Rechargeable battery

- This appliance contains batteries that are non-replaceable.
- The rechargeable battery is a consumable item (not covered by the product warranty).
- The battery is expected to last 300 charge cycles. If the device no longer operates for more than a short period of time even after the battery has been recharged, the battery may have reached the end of its service life.
- For a battery replacement, please contact customer service.

Disposing of the device

- Follow your local government rules for proper disposal of this product.
- This product uses a lithium-ion battery.



Contacts for overseas customers

- For countries and regions outside Japan, please contact the store where you purchased your product for any inquiries. Alternatively, check the YA-MAN GLOBAL Official Website (<https://www.ya-man.com/en/>)

Official website: <https://www.ya-man.com/en/>

Product Specifications

Rated input	5.0 V DC 0.5 A (charging base)
Rated output	5.0 V DC 0.5 A (charging base)
Power consumption	Approx. 1.25 W (when charging)
Battery	Lithium-ion battery
Charging time	Approx. 1.5 to 2.5 hours *May vary depending on the use environment.
Operating time	Approx. 30 minutes (at room temperature of 28°C when fully charged). The actual charging time will differ depending on the remaining battery level, the charging environment, and the AC adapter (commercially available) being used.
Outside dimensions	Approx. W55 x D26 x H141 (mm) (with Wide head)
Weight	Main unit: Approx. 53 g. Wide Head: Approx. 38 g. Spoon Head: Approx. 34 g
Materials	Main unit: ABS resin Wide head and Spoon head: ABS resin, stainless steel Cover common for Wide head and Spoon head: PC resin
Country of origin	China
Accessories	Wide head (attached to the main unit): One Spoon head: One Cover common for Wide head and Spoon head: One Charging base (with USB cord): One User's Manual

Memo

