


How to use RF Booster Pad

Put 1 sheet of RF Booster Pad onto the treatment head and hold it in place with a cotton stopper.

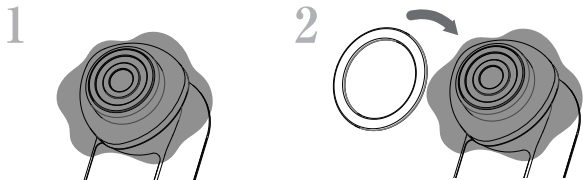

※Always attach or remove the pad while the power is off or in standby mode.

Modes to use with

Deep Cleanse: CLEAN mode

【How to use】 After cleansing your face, wipe off water from your skin before use. After use, do not wipe or wash away. Move onto the next skincare regime.

※After finishing one mode, change to a new sheet if you wish to apply another mode.

Clear Deep Moisture : MOISTURE mode, EMS UP mode, RF mode

【How to use】 Use on clean skin. After use, allow the left-over skincare to blend with your skin.

※You may use the same sheet for two modes without changing to a new one.